



## LumiVibe Go2/Go3 Operations Guide

All Rights Reserved:

No part of this guide covered by the copyright herein may be reproduced, transmitted, stored or used in any form or by any means graphic, electronic, or mechanical, including but not limited to photocopying, recording, scanning, digitizing, taping, web distribution, information networks or information storage and retrieval systems without the prior written permission of LumiCeuticals.

# Chapter 1 Welcome to LumiCeuticals

Welcome to our LumiCeuticals Light Owner Community. We are so excited to have you join us. We see such a wonderful future for harmonic light as the ultimate self-healing tool. We are thrilled to be Lighting The Way....

Lighting the Way for pain-free living with FDA Clearance for Pain Relief and Injury Recovery

Lighting the Way for your power to recover and thrive.

Lighting the Way for your family, pets and community to shine with health and vitality

Lighting the Way for optimum performance and recovery

Lighting the Way for peace and harmony in mind, body and spirit.

Lighting the Way for your inner radiance to shine.

Lighting the Way for a bright financial future for all.

We invite you to use and share your LumiCeuticals light energy system every day.

## 1. The LumiCeuticals Advantage:

- Our systems are all made in an FDA and ISO-9000 approved zero waste medical manufacturing facility and your whole system is backed by a three year full warranty and service in the USA.
- Our systems have medical grade invisium diodes in specific wavelengths of red, infrared and blue light that are well researched for safety and efficacy.
- Our systems combine light and vibrational energy by offering the important added benefit of pulsing the light at 20 highly accurate standard frequencies and a possible 9999 custom frequency options for a profound full body light energy experience that is simple to use and program and unique in the marketplace.
- Each controller has an automatic timer and contains EMF reduction chips for your health and safety.
- All LumiCeuticals Systems come with additional Educational opportunities to further support your use of your harmonic light system. For Further Information Contact:  
[marta@shinewithlight.com](mailto:marta@shinewithlight.com)

# LumiVIBE Go2/3 Operations Guide

Thank you for choosing LumiCeuticals Light equipment. This portable, re-chargeable hand-held light unit will be one of your favorite tools for creating a healthy environment for your body. It is light-weight, easy to use and effective for spot treatments and quick sessions. Please read this guide thoroughly prior to operations in order to maximize the life and efficiency of the product and to provide you with better results. This guide is meant to provide necessary information for: operations, warnings, maintenance and precautionary practices, and warranty regarding this equipment.

**Disclaimer:** LumiCeuticals makes no claims regarding the ability of this or any LumiCeuticals product to cure, mitigate, fix, prevent, repair or heal any disease. Information herein should not be in lieu of the advice of a medical professional. This guide is for educational purposes only. By reading it, you agree to all terms and conditions.

## **Indications for Use:**

Your LumiVibe Go2/Go3 system provides heat and light which studies have shown to relieve pain, stiffness and muscle spasm as well and increasing local blood circulation and relief of injury.

Your LumiVIBE Go2 system will provide many hours of convenient light therapy. This system is ideal for travel or for areas where electrical plugs are not available. The cordless hand held portion of the unit provides a deep, targeted light experience that is highly effective. The blues buster pad and eye mask offer additional light application possibilities. The ability to use variable factory set frequencies makes this a perfect tool for your wellness toolbox. Use it anywhere: in your favorite chair, in your hotel room, car, or campground. We encourage you to take your LumiVIBE Go2 or Go3 everywhere you go!

## **Your System May Contain:**

### **1 LumiVIBE Go2 Portable Handheld Controller**



The 'star' design emits 12 red and 12 near infrared light emitting diodes and can safely be used anywhere on the body except directly into the open eyes. There are 10 light pulsation settings. Settings 1-7 are each individual Nogier harmonic resonance frequencies. Setting A pulses the Nogier settings 1-7 sequentially for 3 minutes each. Setting B sweeps Solfeggio frequencies; Setting C sweeps a Pythagorean 3-2 ratio series of 7 frequencies. See setting guides for details. Do not shine directly into the eyes unless directed to do so by a licensed medical practitioner. Advance the setting by holding down the top button until you see the display advance to your desired setting.


### **1 Controller Carrying Case**

**1 LumiVIBE Go2 Blues Buster Portable Pad:** 50 Near Infrared and 40 Blue Light Diodes. Attaches exclusively to the hand-held controller, can provide light with or without the hand-held controller head illuminated. Not designed for use directly over the eyes.



**1 LumiVibe Go2 Red/Blue Relaxation Eye Mask-** 52 Blue and 52 Red Diodes.  
and/or

**1 LumiVibe Go2 Red/Yellow/IR Restoration Eye Mask** - 27 Red, 44 Near Infrared, 27 Yellow Diodes  
Both attach exclusively to the hand-held controller, can be used over closed eyes and can be used to apply light elsewhere on the body as well providing light with or without the hand-held controller illuminated.

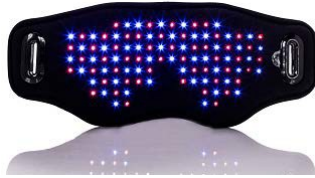

**1 Wall Plug In Charger**

**1 Car or USB Plug In Charger**

Important Note: **The Lumi supplied Go2/Go3 chargers are the only chargers that can be safely used with your controller.**

The portable handheld should arrive fully charged. Should this not be the case or when the charge expires do the following:

**Charging the battery:** Insert either the supplied wall or usb 5V power supply jack into the charging port located on the bottom of the hand-held device. Plug into power source through the car, usb or wall plug in. **Use only the LumiCeuticals provided power cable for this this system.**

- The unit needs charging when the red indicator light on the top right-hand side of the controller is illuminated. A full charge will take 4-6 hours. The green light on the right-hand side of the controller will illuminate when the charge is complete. Note: Attaching a pad to the controller uses more power and will shorten the length of a charge. A charged system will typically have a battery life of 2-3 sessions
- Your unit may be charged in your car. Use the connecting cords provided. LumiVIBE Go2/Go3 is great for travel. Always use with USA certified plugs or use an adaptor.
- LumiVIBE Go2/Go3 may be used while the unit is charging. Take care to avoid electrical shock. Do not leave a charging unit unattended if there are children present.

## **Using Your System:**

**Note: The controller may be used alone or with one pad. The pad may be used with or without the controller illuminated.**

1. **Connecting the pad to the controller:** Visually line up the pins on the blues buster port, before inserting the pad plug into the controller to avoid bending the pins. Be sure that the device is turned off when connecting the Pad. Once the controller is turned on, the pad will light up as well. To use the pad exclusive of the handheld controller lights, depress and hold the channel mode button for 5 seconds to turn off the handheld lights.

2. **Turning on the system:** Press and hold the On/Off button conveniently located mid-way on the controller. Each session will last 20 minutes and will automatically turn off. To turn the unit off at any time, simply press and hold the On/Off button until all lights turn off.

3. **Turning on only the pad:** To turn on the lights on the pad but not the handheld controller, press the top button of the handheld controller and hold it down for 6 seconds. Then advance to the desired session setting by

4. **Changing the settings:** Use the Mode selection button located approximately two inches above the On/Off button. The Mode Display will enable you to know which setting you are choosing. Once the unit is turned on, you can advance to any of the pre-programmed settings. Refer to the setting guides for Nogier Frequency Settings A and 1-7, the Solfeggio Guide for the B setting, and the C setting guide.

## **Caring for Your system:**

### **Do NOT:**

Operate under water

Expose to moisture for a long period of time

Expose to heat and/or sub-freezing temperatures

Crease or excessively bend cords or wrap around pads or device

Stand or walk on unit or apply pressure to any part of the unit.

Stare directly at the light diodes.

### **DO:**

Store in a cool, dry place

Clean gently with a natural, non-alcoholic, non-bleach solution with a soft, damp cloth.

Use only with the power supply provided specifically for this controller by LumiCeuticals.

Using any power source other than the one provide will damage the controller resulting in a high-pitched noise until the battery runs completely down and will void the product warranty.

## **Contraindications:**

Do not use over womb if pregnant, or over the site of an active cancerous growth.

Do not place the red/infrared controller or blue/infrared buster pad over the eyes without medical approval.

Do not use if light sensitive or taking light sensitive drugs such as tetracycline, without medical approval.

Do not use if prone to seizures or epilepsy.

If you have a medical condition, check with your doctor before using this device.

## **Other Precautions and Considerations:**

Do not look directly at the light diodes.

Do not use over topical creams or ointments: they may increase the risk of skin irritation

Not recommended for simultaneous use with ice or cryotherapy.

Monitor physical reaction to light session to prevent irritation or burning sensation.

Use only with building code approved electrical outlets.

Do not use in water, on a wet surface, or when skin is wet.

Do not use pins or metallic materials to hold the light pad in place.

When session is finished, press and hold the 'off' button to turn the power off to the LEDS

Disconnect the power plug gently to prevent damage to the unit.

Cease use if an area of the skin becomes sensitive or too warm.

Apply the light with gentle pressure; pressing harder does not render the lights more effective.

## **Warranty Information:**

LumiCeuticals warrants each piece of the equipment to be free of defects for a period of three years. An optional extended warranty for two additional years is available within the first 30 days of purchase. The warranty does not extend to any liability for medical expenses or for any other direct, indirect, or consequential damages caused by failure, defect or malfunction of any LumiCeuticals equipment.

## **Warranty and Repair:**

Please return defective or non-working equipment to a LumiCeuticals authorized service center. Please notify us right away if you should notice any performance issue by filling out the warranty/service request form at [www.shinewithlight.com/contact us](http://www.shinewithlight.com/contact-us) . LumiCeuticals will contact you within 1 business day after they receive your customer service report.